

聯合骨科器材股份有限公司

會議記錄單

頁次: 1 / 3

時 間：2018 年 04 月 03 日 (星 期 二) 13 時 00 分 ~ 13 時 30 分					
地 點：HC_55A, TP_A					
主 席：倪仲達 記 錄：倪仲達					
會議名稱：Conformity stem 原型確認會議 權責單位：UR3					
出席人員：劉育良、陳文全、陸正光、鄭雅文、倪仲達、簡惠蓮、何浩君、嶺榮娟、羅正展、林師緯、劉雅甄、張宗維、鄭裕耀、鍾宇琪、林成遠、邱鈺芬					
項次	會 議 內 容	負 責 人	預計日期	實際完成日	追蹤結果
1	UPD建議cemented 與cementless的CE應一併送件，請評估cemented stem的開發與送件時間	倪仲達	4/20		
2	與UR4法規確認 TFDA工單樣品的數量為10支	倪仲達	4/13	4/13	
3	產品預計申請地區為： CE , FDA, TFDA, CFDA ，請UPD協助確認是否有日本送件需求	UPD			
4	實質對等產品為Depuy 的Corail stem	--			
5	確認所有輸出文件的時間， 經確認後所有輸出文件時間為：		倪仲達	4/11	
	Stage 3 設計輸出	完成日期			
	專案計劃書	2018/1/2			
	臨床評估計畫書/報告書	2018/2/14			
	醫療器材基本安全要求檢查表	2018/1/15			
	風險管理計劃書	2018/3/26			
	風險管理報告	2018/3/26			
	產品設計輸入/產出檢查表	2018/4/3			
	產品系統明細表	2018/3/2			
產品測試表	2018/4/2				

聯合骨科器材股份有限公司

會議記錄單

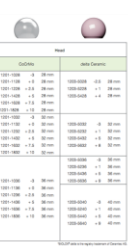
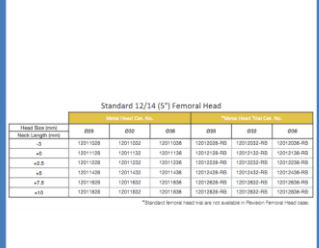
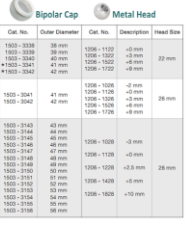
頁次: 2 / 3

6	修正產品名，除Conformity 大寫外，其餘小寫 修正後如下:					倪仲達	4/11		
	品號	品名	中文品名	Stem length	Offset	Distal width			
	11101001	Conformity stem, collared, #1	優康有領股骨柄, #1	115.0 mm	36.0mm	7.5mm			
	11101002	Conformity stem, collared, #2	優康有領股骨柄, #2	119.5 mm	36.5 mm	8.1 mm			
	11101003	Conformity stem, collared, #3	優康有領股骨柄, #3	124.0 mm	37.5 mm	8.5 mm			
	11101004	Conformity stem, collared, #4	優康有領股骨柄, #4	128.5 mm	38.0 mm	9.5 mm			
	11101005	Conformity stem, collared, #5	優康有領股骨柄, #5	133.0 mm	39.0 mm	10.0 mm			
	11101006	Conformity stem, collared, #6	優康有領股骨柄, #6	137.5 mm	39.5mm	10.0 mm			
	11101007	Conformity stem, collared, #7	優康有領股骨柄, #7	142.0 mm	40.0mm	10.0 mm			
	11101008	Conformity stem, collared, #8	優康有領股骨柄, #8	146.5 mm	41.0 mm	10.0 mm			
	11101009	Conformity stem, collared, #9	優康有領股骨柄, #9	151.0 mm	41.5 mm	10.0 mm			
	11101010	Conformity stem, collared, #10	優康有領股骨柄, #10	155.5 mm	42.5 mm	11.0 mm			
	11101011	Conformity stem, collared, #11	優康有領股骨柄, #11	160.0 mm	43.5 mm	11.0 mm			
	11101201	Conformity stem, high offset, collared, #1	優康有領高偏距股骨柄, #1	115.0 mm	43.0 mm	7.5mm			
	11101202	Conformity stem, high offset, collared, #2	優康有領高偏距股骨柄, #2	119.5 mm	43.5 mm	8.1 mm			
	11101203	Conformity stem, high offset, collared, #3	優康有領高偏距股骨柄, #3	124.0 mm	44.5mm	8.5 mm			
	11101204	Conformity stem, high offset, collared, #4	優康有領高偏距股骨柄, #4	128.5 mm	45.0 mm	9.5 mm			
	11101205	Conformity stem, high offset, collared, #5	優康有領高偏距股骨柄, #5	133.0 mm	46.0 mm	10.0 mm			
	11101206	Conformity stem, high offset, collared, #6	優康有領高偏距股骨柄, #6	137.5 mm	46.5 mm	10.0 mm			
	11101207	Conformity stem, high offset, collared, #7	優康有領高偏距股骨柄, #7	142.0 mm	47.0 mm	10.0 mm			
	11101208	Conformity stem, high offset, collared, #8	優康有領高偏距股骨柄, #8	146.5 mm	48.0mm	10.0 mm			
	11101209	Conformity stem, high offset, collared, #9	優康有領高偏距股骨柄, #9	151.0 mm	48.5 mm	10.0 mm			
	11101210	Conformity stem, high offset, collared, #10	優康有領高偏距股骨柄, #10	155.5 mm	49.5 mm	11.0 mm			
	11101211	Conformity stem, high offset, collared, #11	優康有領高偏距股骨柄, #11	160.0 mm	50.5 mm	11.0 mm			
	品號	品名	中文品名	Stem length	Offset	Distal width			
	11103001	Conformity stem, #1	優康股骨柄, #1	115.0 mm	36.0mm	7.5mm			
	11103002	Conformity stem, #2	優康股骨柄, #2	119.5 mm	36.5 mm	8.1 mm			
	11103003	Conformity stem, #3	優康股骨柄, #3	124.0 mm	37.5 mm	8.5 mm			
	11103004	Conformity stem, #4	優康股骨柄, #4	128.5 mm	38.0 mm	9.5 mm			
	11103005	Conformity stem, #5	優康股骨柄, #5	133.0 mm	39.0 mm	10.0 mm			
	11103006	Conformity stem, #6	優康股骨柄, #6	137.5 mm	39.5mm	10.0 mm			
	11103007	Conformity stem, #7	優康股骨柄, #7	142.0 mm	40.0mm	10.0 mm			
	11103008	Conformity stem, #8	優康股骨柄, #8	146.5 mm	41.0 mm	10.0 mm			
	11103009	Conformity stem, #9	優康股骨柄, #9	151.0 mm	41.5 mm	10.0 mm			
	11103010	Conformity stem, #10	優康股骨柄, #10	155.5 mm	42.5 mm	11.0 mm			
	11103011	Conformity stem, #11	優康股骨柄, #11	160.0 mm	43.5 mm	11.0 mm			
	11103201	Conformity stem, high offset, #1	優康高偏距股骨柄, #1	115.0 mm	43.0 mm	7.5mm			
	11103202	Conformity stem, high offset, #2	優康高偏距股骨柄, #2	119.5 mm	43.5 mm	8.1 mm			
	11103203	Conformity stem, high offset, #3	優康高偏距股骨柄, #3	124.0 mm	44.5mm	8.5 mm			
	11103204	Conformity stem, high offset, #4	優康高偏距股骨柄, #4	128.5 mm	45.0 mm	9.5 mm			
	11103205	Conformity stem, high offset, #5	優康高偏距股骨柄, #5	133.0 mm	46.0 mm	10.0 mm			
	11103206	Conformity stem, high offset, #6	優康高偏距股骨柄, #6	137.5 mm	46.5 mm	10.0 mm			
	11103207	Conformity stem, high offset, #7	優康高偏距股骨柄, #7	142.0 mm	47.0 mm	10.0 mm			
	11103208	Conformity stem, high offset, #8	優康高偏距股骨柄, #8	146.5 mm	48.0mm	10.0 mm			
	11103209	Conformity stem, high offset, #9	優康高偏距股骨柄, #9	151.0 mm	48.5 mm	10.0 mm			
	11103210	Conformity stem, high offset, #10	優康高偏距股骨柄, #10	155.5 mm	49.5 mm	11.0 mm			
	11103211	Conformity stem, high offset, #11	優康高偏距股骨柄, #11	160.0 mm	50.5 mm	11.0 mm			
	11105202	Conformity stem, coxa vara, collared, #2	優康有領內翻股骨柄, #2	119.5 mm	43.5 mm	8.1 mm			
	11105203	Conformity stem, coxa vara, collared, #3	優康有領內翻股骨柄, #3	124.0 mm	44.5mm	8.5 mm			
	11105204	Conformity stem, coxa vara, collared, #4	優康有領內翻股骨柄, #4	128.5 mm	45.0 mm	9.5 mm			
	11105205	Conformity stem, coxa vara, collared, #5	優康有領內翻股骨柄, #5	133.0 mm	46.0 mm	10.0 mm			
	11105206	Conformity stem, coxa vara, collared, #6	優康有領內翻股骨柄, #6	137.5 mm	46.5 mm	10.0 mm			
	11105207	Conformity stem, coxa vara, collared, #7	優康有領內翻股骨柄, #7	142.0 mm	47.0 mm	10.0 mm			
	11105208	Conformity stem, coxa vara, collared, #8	優康有領內翻股骨柄, #8	146.5 mm	48.0mm	10.0 mm			
	11105209	Conformity stem, coxa vara, collared, #9	優康有領內翻股骨柄, #9	151.0 mm	48.5 mm	10.0 mm			
	11105210	Conformity stem, coxa vara, collared, #10	優康有領內翻股骨柄, #10	155.5 mm	49.5 mm	11.0 mm			
	11105211	Conformity stem, coxa vara, collared, #11	優康有領內翻股骨柄, #11	160.0 mm	50.5 mm	11.0 mm			

聯合骨科器材股份有限公司

會議記錄單

頁次: 3 / 3

7	搭配產品： 預計搭配球頭 Ceramic head  Metal head  Bipolar 	--																														
7	測試項目及允收標準： <table><tr><th>測試項目</th><th>樣品規格</th><th>允收標準</th></tr><tr><td>Range of Motion analysis</td><td> •Acetabular cup liner: PE 1403-1044 (XPE cup liner, ID 28mm, 10°x44mm) •Femoral Head: PE 1203-5428 Delta ceramic femoral heads. •Stem: HO#1 </td><td> ISO 21535 (100° in flexion-extension, 90° in internal/external rotation, 60° in abduction/adduction) </td></tr><tr><td>Stem fatigue analysis and test</td><td>•The Worst case of ANSYS simulations</td><td>ISO 7206-4 Pass</td></tr><tr><td>Neck fatigue analysis and test</td><td> •The Worst case of ANSYS simulations •Stem: conformity stem high offset •Femoral head, Ø 36mm, +10mm (PE 12011836) </td><td>ISO 7206-6 Pass</td></tr><tr><td>Pull-off Strength Testing</td><td>•引用UOC-RD-TR-11015</td><td>ISO 7206-10</td></tr><tr><td>Impingement analysis and test</td><td>•The Worst case of ANSYS simulations</td><td> ASTM F2582 not fracture when subjected to 1 million cycles impingement test. </td></tr><tr><td>內毒素</td><td> •The highest surface area sample •CHO #11 </td><td>ISO 11137</td></tr><tr><td>滅菌確效</td><td>•滅菌確效評估表</td><td>ISO11737-1</td></tr><tr><td>HA Coating Characterization Testing</td><td>MAF of coating company (Euro-coating)</td><td>See attachment.</td></tr></table>	測試項目	樣品規格	允收標準	Range of Motion analysis	•Acetabular cup liner: PE 1403-1044 (XPE cup liner, ID 28mm, 10°x44mm) •Femoral Head: PE 1203-5428 Delta ceramic femoral heads. •Stem: HO#1	ISO 21535 (100° in flexion-extension, 90° in internal/external rotation, 60° in abduction/adduction)	Stem fatigue analysis and test	•The Worst case of ANSYS simulations	ISO 7206-4 Pass	Neck fatigue analysis and test	•The Worst case of ANSYS simulations •Stem: conformity stem high offset •Femoral head, Ø 36mm, +10mm (PE 12011836)	ISO 7206-6 Pass	Pull-off Strength Testing	•引用UOC-RD-TR-11015	ISO 7206-10	Impingement analysis and test	•The Worst case of ANSYS simulations	ASTM F2582 not fracture when subjected to 1 million cycles impingement test.	內毒素	•The highest surface area sample •CHO #11	ISO 11137	滅菌確效	•滅菌確效評估表	ISO11737-1	HA Coating Characterization Testing	MAF of coating company (Euro-coating)	See attachment.	--			
測試項目	樣品規格	允收標準																														
Range of Motion analysis	•Acetabular cup liner: PE 1403-1044 (XPE cup liner, ID 28mm, 10°x44mm) •Femoral Head: PE 1203-5428 Delta ceramic femoral heads. •Stem: HO#1	ISO 21535 (100° in flexion-extension, 90° in internal/external rotation, 60° in abduction/adduction)																														
Stem fatigue analysis and test	•The Worst case of ANSYS simulations	ISO 7206-4 Pass																														
Neck fatigue analysis and test	•The Worst case of ANSYS simulations •Stem: conformity stem high offset •Femoral head, Ø 36mm, +10mm (PE 12011836)	ISO 7206-6 Pass																														
Pull-off Strength Testing	•引用UOC-RD-TR-11015	ISO 7206-10																														
Impingement analysis and test	•The Worst case of ANSYS simulations	ASTM F2582 not fracture when subjected to 1 million cycles impingement test.																														
內毒素	•The highest surface area sample •CHO #11	ISO 11137																														
滅菌確效	•滅菌確效評估表	ISO11737-1																														
HA Coating Characterization Testing	MAF of coating company (Euro-coating)	See attachment.																														
8	導量產數量： <table><tr><th>測試</th><th>測試樣品</th></tr><tr><td>Stem fatigue</td><td>1110-3201(6 pieces)</td></tr><tr><td>Neck fatigue</td><td>1110-3201(6 pieces)</td></tr><tr><td>內毒素</td><td>1110-1211(6 pieces)</td></tr><tr><td>滅菌確效</td><td>1110-1211(10 pieces)</td></tr></table>	測試	測試樣品	Stem fatigue	1110-3201(6 pieces)	Neck fatigue	1110-3201(6 pieces)	內毒素	1110-1211(6 pieces)	滅菌確效	1110-1211(10 pieces)	--																				
測試	測試樣品																															
Stem fatigue	1110-3201(6 pieces)																															
Neck fatigue	1110-3201(6 pieces)																															
內毒素	1110-1211(6 pieces)																															
滅菌確效	1110-1211(10 pieces)																															
9	提供標籤狀況給UR4確認，詳如附件。	倪仲達	4/16																													
10	包裝部分等競品購入之後再進行確認	倪仲達	4/16																													
11	Stem fatigue 的部分以distal width 為7.5mm進行測試	--																														